Mitch TRH Total Conventional Hip Investigation

Note: This analysis compares the Mitch TRH acetabular prosthesis with all other total conventional hip prostheses.

This prosthesis has been identified as having a significantly higher rate of revision. For a detailed explanation of the process used by the Registry that results in identification of prostheses that have a higher than anticipated rate of revision please refer to the Prostheses with Higher than Anticipated Rates of Revision chapter of the most recent AOANJRR Annual Report, https://aoanjrr.sahmri.com/annual-reports-2024.

Note: Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator. Procedures using prostheses with no recorded use in 2023 are excluded from the comparator.

TABLE 1

Revision Rate of Primary Total Conventional Hip Replacement

The revision rate of the Mitch TRH total conventional hip prosthesis is compared to all other total conventional hip prostheses.

Table 1: Revision Rates of Primary Total Conventional Hip Replacement

Component	N Revised	N Total	Obs. Years	Revisions/100 Obs. Yrs (95% Cl)
Mitch TRH	176	731	8604	2.05 (1.75, 2.37)
Other Total Conventional Hip	19249	538541	3454641	0.56 (0.55, 0.57)
TOTAL	19425	539272	3463245	0.56 (0.55, 0.57)

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Mitch TRH total conventional hip prosthesis is compared to all other total conventional hip prostheses.

CPR	1 Yr	2 Yrs	3 Yrs	4 Yrs	5 Yrs	6 Yrs	7 Yrs	8 Yrs
Mitch TRH	1.5 (0.8, 2.7)	3.2 (2.1, 4.7)	4.6 (3.3, 6.4)	6.2 (4.6, 8.2)	7.7 (6.0, 10.0)	9.4 (7.4, 11.8)	10.9 (8.8, 13.4)	12.4 (10.1, 15.1)
Other Total Conventional Hip	1.7 (1.7, 1.8)	2.2 (2.1, 2.2)	2.5 (2.5, 2.5)	2.8 (2.7, 2.8)	3.1 (3.0, 3.1)	3.3 (3.3, 3.4)	3.6 (3.6, 3.7)	3.9 (3.8, 4.0)
CPR	9 Yrs	10 Yrs	11 Yrs	12 Yrs	13 Yrs	14 Yrs	15 Yrs	16 Yrs
Mitch TRH	13.5 (11.1, 16.3)	15.6 (13.1, 18.6)	17.9 (15.2, 21.1)	21.4 (18.4, 24.7)	23.0 (19.9, 26.4)	25.5 (22.2, 29.2)	27.3 (23.8, 31.1)	29.8 (26.0, 34.0)
Other Total Conventional Hip	4.2 (4.2, 4.3)	4.5 (4.5, 4.6)	4.9 (4.8, 5.0)	5.3 (5.2, 5.4)	5.7 (5.6, 5.8)	6.1 (5.9, 6.2)	6.5 (6.3, 6.6)	6.9 (6.7, 7.0)
CPR	17 Yrs	18 Yrs	19 Yrs	20 Y	′rs 21	Yrs	22 Yrs	23 Yrs
Mitch TRH								
Other Total Conventional Hip	7.3 (7.1, 7.4)) 7.7 (7.5, 7.5	8) 8.2 (8.0,	8.4) 8.5 (8.	2, 8.7) 9.0 (8.7, 9.3) 9.7	(9.2, 10.1)	10.3 (9.5, 11.2)


Table 2: Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

FIGURE 1

Yearly Cumulative Percent Revision of Primary Total Conventional Hip Replacement

The yearly cumulative percent revision of the Mitch TRH total conventional hip prosthesis is compared to all other total conventional hip prostheses. In addition, hazard ratios are reported.

Hazard ratios are reported for specific time periods during which the hazard ratio is constant. This is done to enable more specific and valid comparisons of the risk of revision over time. The pattern of variation in risk has important implications with respect to the underlying reasons for any difference.

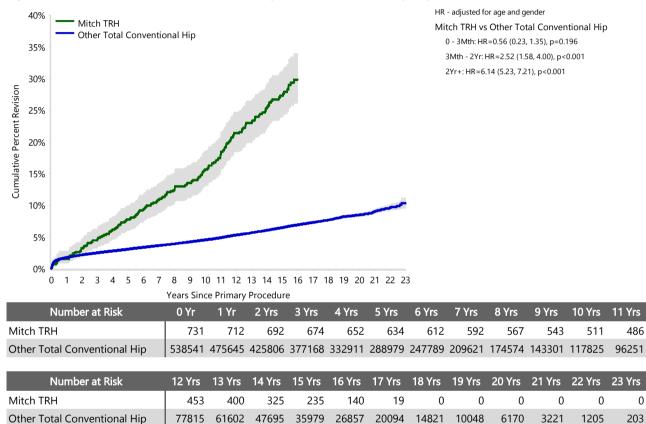



Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

This table identifies the diagnosis of the primary procedure which was subsequently revised. This information is provided as there is a variation on outcome depending on the primary diagnosis. It is therefore important when considering the reasons for a higher than anticipated rate of revision that there is identification of the primary diagnosis. This information should be compared to the primary diagnosis for the revisions of all other total conventional hip prostheses.

Table 3: Primary Diagnosis for Revised Primary Total Conventional Hip Replacement

	Mitch	n TRH	Other Total Co	nventional Hip
Primary Diagnosis	Number	Percent	Number	Percent
Osteoarthritis	159	90.3	15951	82.9
Fractured Neck Of Femur	4	2.3	1420	7.4
Osteonecrosis	9	5.1	859	4.5
Developmental Dysplasia			320	1.7
Rheumatoid Arthritis	1	0.6	208	1.1
Failed Internal Fixation	1	0.6	151	0.8
Tumour			149	0.8
Other Inflammatory Arthritis	2	1.1	106	0.6
Fracture/Dislocation			53	0.3
Other			17	0.1
Arthrodesis Takedown			15	0.1
TOTAL	176	100.0	19249	100.0

Reasons for Revision

This is reported in two ways: a percentage of primary procedures revised and as a percentage of all revision procedures.

% Primaries Revised: This shows the proportional contribution of each revision diagnosis as a percentage of the total number of primary procedures. This percentage can be used to approximate the risk of being revised for that diagnosis. Differing percentages between groups, with the same distribution of follow up time, may identify problems of concern.

% Revisions: The number of revisions for each diagnosis is expressed as a percentage of the total number of revisions. This shows the distribution of reasons for revision within a group but cannot be used as a comparison between groups.

Table 4: Primary T	Total Conventional Hip	Replacement -	Reason for Revision	(Follow-u	p Limited to 17.5 Years)
rable in rinnary i	otal controllational inp	replacement		(1011011 01	

		Mitch TRH		Othe	r Total Convention	al Hip
Revision Diagnosis	Number	% Primaries Revised	% Revisions	Number	% Primaries Revised	% Revisions
Infection	25	3.4	14.2	4486	0.8	23.5
Prosthesis Dislocation/Instability	3	0.4	1.7	4362	0.8	22.9
Fracture	9	1.2	5.1	4237	0.8	22.2
Loosening	35	4.8	19.9	3774	0.7	19.8
Pain	4	0.5	2.3	335	0.1	1.8
Leg Length Discrepancy	2	0.3	1.1	291	0.1	1.5
Malposition	1	0.1	0.6	266	0.0	1.4
Lysis	15	2.1	8.5	202	0.0	1.1
Implant Breakage Stem	23	3.1	13.1	188	0.0	1.0
Implant Breakage Acetabular Insert				130	0.0	0.7
Incorrect Sizing				103	0.0	0.5
Wear Acetabular Insert				102	0.0	0.5
Metal Related Pathology	58	7.9	33.0	77	0.0	0.4
Implant Breakage Acetabular				71	0.0	0.4
Wear Head				46	0.0	0.2
Tumour				44	0.0	0.2
Implant Breakage Head				32	0.0	0.2
Heterotopic Bone				26	0.0	0.1
Wear Acetabulum				10	0.0	0.1
Osteonecrosis				2	0.0	0.0
Progression Of Disease				2	0.0	0.0
Synovitis				1	0.0	0.0
Other	1	0.1	0.6	292	0.1	1.5
N Revision	176	24.1	100.0	19079	3.5	100.0
N Primary	731			538541		

Note: This table is restricted to revisions within 17.5 years for all groups to allow a time-matched comparison of revisions. Note: Prostheses no longer used in 2023 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

FIGURE 2

Cumulative Incidence Revision Diagnosis of Primary Total Conventional Hip Replacement

This figure details the cumulative incidence of the most common reasons for revision. The five most common reasons for revision are included as long as each of these reasons account for more than 10 procedures or at least 5% of all revisions for the Mitch TRH total conventional hip prosthesis. A comparative graph is provided of the cumulative incidence for the same reasons for revisions for all other total conventional hip prostheses.

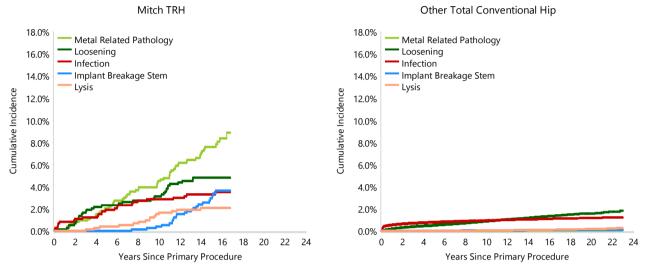


Figure 2: Cumulative Incidence Revision Diagnosis for Primary Total Conventional Hip Replacement

Type of Revision Performed for Primary Total Conventional Hip Replacement

This analysis identifies the components used in the revision of the Mitch TRH total conventional hip prosthesis and compares it to the components used in the revision of all other total conventional hip prostheses.

The reason this analysis is undertaken is to identify whether there is one or more components which are being replaced that differ from the components replaced for revisions of all other total conventional hip prostheses i.e. is there a difference in the type of revision undertaken for the Mitch TRH total conventional hip prosthesis compared to all other total conventional hip prostheses.

	Mitc	h TRH	Other Total Co	nventional Hip
Type of Revision	Number	Percent	Number	Percent
Femoral Component	30	17.0	6313	33.1
Acetabular Component	45	25.6	3559	18.7
THR (Femoral/Acetabular)	81	46.0	2214	11.6
Cement Spacer	13	7.4	620	3.2
Removal of Prostheses			98	0.5
Reinsertion of Components			28	0.1
Total Femoral			9	0.0
Bipolar Head and Femoral			7	0.0
Saddle			1	0.0
N Major	169	96.0	12849	67.3
Head/Insert	5	2.8	4806	25.2
Head Only	2	1.1	925	4.8
Vinor Components			309	1.6
nsert Only			186	1.0
Bipolar Only			2	0.0
Cement Only			1	0.0
Head/Neck			1	0.0
N Minor	7	4.0	6230	32.7
TOTAL	176	100.0	19079	100.0

Note: This table is restricted to revisions within 17.5 years for all groups to allow a time-matched comparison of revisions. Note: Prostheses no longer used in 2023 are excluded from the comparator. Procedures using metal/metal prostheses with head size larger than 32mm are excluded from the comparator.

Revision Rates of Mitch TRH Primary Total Conventional Hip Replacement by Fixation

This analysis is provided as some prostheses have more than one fixation option. Additionally there are prostheses where an alternative to the recommended approach to fixation was used e.g. a cementless prosthesis that has been cemented or vice-versa.

Table 6: Revised Number of Mitch TRH Primary Total Conventional Hip Replacement by Fixation

Fixation	N Revised	N Total
Cementless	128	429
Hybrid (Femur Cemented)	48	302
TOTAL	176	731

TABLE 7

Revision Rates of Mitch TRH Primary Total Conventional Hip Replacement by Bearing Surface

This analysis is provided as some prostheses are combined with a variety of bearing surfaces. All bearing surfaces used with this prosthesis are listed.

Table 7: Revised Number of Mitch TRH Primary Total Conventional Hip Replacement by Bearing Surface

Bearing Surface	N Revised	N Total
Metal/Metal	176	730
Unknown	0	1
TOTAL	176	731

Revision Rates of Primary Total Conventional Hip Replacement by State

This enables a state by state variation to be identified for the Mitch TRH total conventional hip prosthesis and provides the comparative data for each of the states for all other total conventional hip prostheses.

The purpose of this analysis is to determine if the higher than anticipated rate of revision has widespread distribution between states. If there is widespread distribution then the reason for the higher than anticipated rate of revision is unlikely to be surgeon specific. If the prosthesis has been used in only a small number of states it is not possible to distinguish if the higher than anticipated rate of revision is related to the prosthesis, surgeon, technique or patient.

Component	State	N Revised	N Total	
Mitch TRH	NSW	70	244	
	VIC	7	20	
	QLD	18	121	
	WA	53	184	
	SA	8	100	
	TAS	5	13	
	ACT/NT	15	49	
Other Total Conventional Hip	NSW	5207	157707	
	VIC	4859	140681	
	QLD	3791	94582	
	WA	2546	63265	
	SA	1825	49398	
	TAS	452	18198	
	ACT/NT	569	14710	
TOTAL		19425	539272	

Table 8: Revised Number of Primary Total Conventional Hip Replacement by State

Number of Revisions of Mitch TRH Primary Total Conventional Hip Replacement by Year of Implant

This analysis details the number of prostheses reported each year to the Registry for the Mitch TRH total conventional hip prosthesis. It also provides the subsequent number of revisions of the primaries reported in that year.

Primary procedures performed in later years have had less follow up time therefore the number revised is expected to be less than the number revised in earlier years. For example, a primary procedure performed in 2023 has a maximum of one year to be revised, whereas a primary procedure performed in 2021 has a maximum of three years to be revised.

Table 9: Number of Revisions of Mitch TRH Primary Total Conventional Hip Replacement by Year of Implant

Year of Implant	Number Revised	Total Number
2006	12	45
2007	70	273
2008	38	164
2009	32	130
2010	22	82
2011	2	37
TOTAL	176	731

Revision Rates of Mitch TRH Primary Total Conventional Hip Replacement by Catalogue Number Range

Many prostheses have a number of catalogue ranges. The catalogue range is specific to particular design features; more than one catalogue range usually indicates a minor difference in design in a particular Mitch TRH prosthesis.

This analysis has been undertaken to determine if the revision rate varies according to the catalogue number range.

Model	Catalogue Range	Catalogue Description	Cement	Material	Coating
Acetabular					
Mitch TRH	MAC99883844-MAC99885864	TITANIUM HA STANDARD RS ACETABULAR CUP	NO	METAL	HA COATED

Table 10: Revised Number of Mitch TRH Primary Total Conventional Hip Replacement by Catalogue Number Range

Acetabular Range	N Revised	N Total
MAC99883844-MAC99885864	176	731
TOTAL	176	731

Revision Rates of Mitch TRH Primary Total Conventional Hip Replacement by Component

A prosthesis may be combined with multiple components. This analysis has been undertaken to determine if the revision rate varies according to the component with which it is combined.

Table 11: Revised Number of Mitch TRH Primary Total Conventional Hip Replacement by Femoral Stem Component

Femoral Stem Component	N Revised	N Total
ABGII	7	21
ABGII (exch neck)	1	2
Accolade I	117	401
Citation	2	4
Exeter V40	47	300
Restoration	2	3
TOTAL	176	731